

NUTRICIA
flocare *Infinity*[™]

REDESIGNED FLOCARE STAND



NOW
MORE FLEXIBLE
THAN EVER

Danone Specialized Nutrition
Danone Place Amsterdam
Taurusavenue 167
2132 LS Hoofddorp, The Netherlands

For more information on your Flocare pump and access to the complete portfolio of feeding tubes, sets and pumps, as well as accessories, training videos, safety information and instructions for use, please visit the Flocare website: www.nutriciaflocare.com

Copyright © 2022 Danone Specialized Nutrition. This material/document 400048_A could include inaccuracies and errors.

NUTRICIA
flocare *Infinity*[™]

REDESIGNED FLOCARE STAND NOW MORE FLEXIBLE THAN EVER

The improved design of the Flocare Infinity Universal Frame means that the Infinity III as well as Infinity II pump can be used anywhere, from hospitals and schools to patient's homes. It can move from the bedside to the tabletop, or further afield in a bag or backpack. It allows users to be fully mobile and independent.

ADVANCED, REDESIGNED STAND

- Improved design integrates the Flocare Infinity stand and the Flocare GO frame into a single, universal frame
- Multifunctional: suitable for tabletop use, use as a bedstand, or use as a mobile frame in a dedicated bag or backpack
- Usable with giving sets equipped with and without dripchamber
- The Universal Frame Adapter allows for quick access to the pump to use it as handheld, tabletop reader or to secure it safely into the frame



TABLETOP/BEDSTAND:

There are 3 'high shelf' positions to facilitate the connection between various reservoirs and giving sets WITH dripchamber.



TABLETOP/BEDSTAND:

2 Lowest shelf positions ensure usability for giving sets without a dripchamber




BACKPACK:

Giving set without dripchamber

TRUSTED FEATURES YOU WILL RECOGNIZE FROM THE INFINITY GO FRAME

- Four rubber foot pads enhance stability and reduce vibration and noise
- Closed base collects spills
- The new frame fits with the existing Flocare backpack range
- Robust frame secures the Flocare Infinity pumps (II and III) and the feeding reservoir while providing easy access to all of its functionalities

NEW AND IMPROVED FEATURES

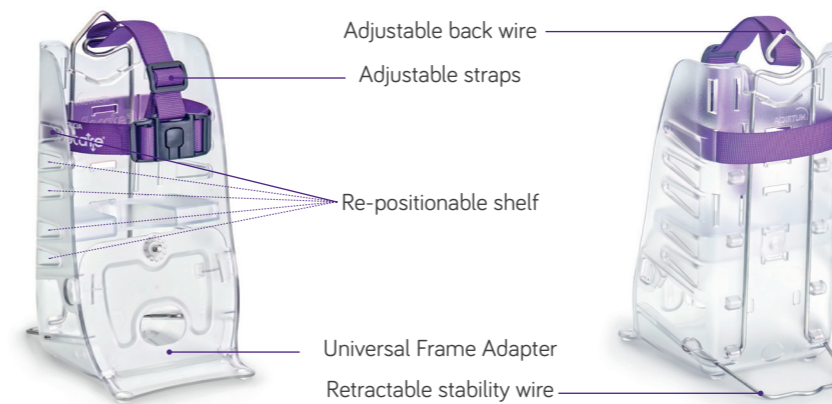

ADJUSTABLE VERTICAL BACK WIRE

- Robust and simple to use design stabilising large reservoirs
- Height adjustment with easy-to-use snap lock
- Allows frame to be carried safely
- Shaped to allow hanging on hook or infusion stand
- Safe to handle – no sharp corners or points

ADJUSTABLE STRAPS WITH IMPROVED BUCKLE DESIGN

- To accommodate and secure reservoirs of various sizes

RE-POSITIONABLE SHELF

- Five positions allow use of a broader range of giving sets, with and without drip chamber
- Compatible with full range of Nutricia feeding reservoirs including other bottles
- Symmetric shelf design enables use upside down without problems

FRAME ADAPTER

- The Universal Frame Adapter allows for quick access to the pump to use it as handheld, tabletop reader or to secure it safely into the frame
- Fits with a range of pumps (a.o. INF II and INF III)
- Offers a screw nut to fix pumps to the adapter
- Provides a click system to easily click in and out the adapter from the frame
- The frameholder has a tubing guide to prevent tube getting stuck between the adapter and the frame and prevents kinking and pinching of the tube

RETRACTABLE HORIZONTAL STABILITY WIRE

- Re-positionable for extra tabletop stability

Manufactured from durable, reliable and safe material.
Dishwasher safe to 55°C.



Backpack with direct access to the pump via the front zipper



Highest shelf position
Pack giving set with dripchamber



Highest shelf position
Pack & Bottle giving set with dripchamber